

REF 14684 1 X 5 ML BOTTLE, CALIBRATOR

PROPOXYPHENE (PPX) CUTOFF CALIBRATOR

The bottle contains 5.0 mL of human urine based ready-to-use liquid spiked with a known concentration of propoxyphene.

INTENDED USE

The Medica Propoxyphene Cutoff Calibrator is intended for in-vitro diagnostic use for the calibration of the propoxyphene assay, which is used on the EasyRA for the qualitative determination of propoxyphene in urine.

SAFETY PRECAUTIONS AND WARNINGS

1. The human based material from which this product has been derived has not been tested for the presence of antibodies to HIV, the hepatitis antigens, or anti-hepatitis antibodies. This material and all patient samples should be handled and disposed of as though they are capable of transmitting infectious diseases.
2. The Propoxyphene Cutoff Calibrator contains sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azides. When disposing, always flush with copious volumes of water to prevent azide build-up.
3. Do not use the calibrator beyond the expiration date.
4. The bottle containing this calibrator should be tightly capped and refrigerated when not in use to prevent evaporation.

INSTRUCTIONS FOR HANDLING, STORAGE AND STABILITY

The calibrator material is stable until the expiration date on the label if stored at 2° – 8°C. The calibrator should not be frozen or exposed to temperatures above 32°C.

INSTRUCTIONS FOR USE

The calibrator material is ready-to-use. Remove bottle cap. Dispense three drops of material into a sample cup and proceed with instructions in the EasyRA Operator's Manual.

VALUE ASSIGNMENT


The table below contains the assigned values for the Propoxyphene Cutoff Calibrator:

Lot No. XXXXX Expires: XXXX-XX-XX

Analyte:	Assigned Value	Units
Propoxyphene	300	ng/mL



Medica Corporation, 5 Oak Park Drive
Bedford, MA 01730-1413 USA

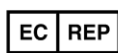

Emergo Europe, Prinsessegracht 20
2514 AP The Hague, The Netherlands

008979-901 R2